

Human GM-CSF, cct-premium

Human Granulocyte Macrophage Colony Stimulating Factor, recombinant

Synonyms: CSF2, GMCSF

PLEASE NOTE: ALWAYS CENTRIFUGE VIAL BEFORE OPENING

Size	Order #	Lot #	Expiry Date
25 µg	1450.950.025		
50 µg	1450.950.050		
100 µg	1450.950.100		
200 µg	1450.950.200		
500 µg	1450.950.500		
1 mg	1450.950.199		

Please enquire for bulk quantities and other vial sizes

Description

Recombinant human Granulocyte Macrophage Colony Stimulating Factor (GM-CSF), a 14,5 kDa protein consisting of 127 amino acid residues (Ala18-Glu144), is a potent species specific stimulator of bone marrow cells and several other cell types. GM-CSF was initially characterized as a growth factor that can support the in vitro colony formation of granulocyte-macrophage progenitors. It is produced by a number of different cell types (including activated T cells, B cells, macrophages, mast cells, endothelial cells and fibroblasts) in response to cytokine or immune and inflammatory stimuli. Besides granulocyte-macrophage progenitors, GM-CSF is also a growth factor for erythroid, megakaryocyte and eosinophil progenitors. On mature hematopoietic cells, GM-CSF is a survival factor for and activates the effector functions of granulocytes, monocytes/macrophages and eosinophils. GM-CSF has also been reported to have a functional role on non-hematopoietic cells. It can induce human endothelial cells to migrate and proliferate. Additionally, GM-CSF can also stimulate the proliferation of a number of tumor cell lines, including osteogenic sarcoma, carcinoma and adenocarcinoma cell lines. GM-CSF is species specific and human GM-CSF has no biological effects on mouse cells. GM-CSF exerts its biological effects through binding to specific cell surface receptors. The high affinity receptors required for human GM-CSF signal transduction have been shown to be heterodimers consisting of a GM-CSF-specific α chain and a common β chain that is shared by the high-affinity receptors for IL-3 and IL-5.

- **Biological Activity** $\geq 1 \times 10^7$ IU/mg
- **Source** *E. Coli*
- **Purity** ≥ 98 % (SDS-PAGE, silver stained)
- **Endotoxin level** ≤ 0.01 ng/µg of protein (≤ 0.1 EU/µg)
- **Stabilizer** None
- **Buffer** PBS, pH 7.2*
- **Physical state** Sterile filtered, lyophilized

Biological Activity

The ED₅₀ of ≤ 0.1 ng/ml was determined by the dose-dependent stimulation of the proliferation of human TF-1 cells. This corresponds to a specific activity of $\geq 1.0 \times 10^7$ units/mg.

Reconstitution

We recommend a quick spin followed by reconstitution in sterile filtered water to a concentration of 0.1 mg/ml. Do not vortex. This solution can then be diluted into other aqueous buffers and stored at 4°C for up to 1 week or at -20°C for future use.

Stability

The lyophilized protein is stable at room temperature for up to 3 weeks and below -20°C until expiry date. Reconstituted human GM-CSF should be stored in working aliquots at -20°C. **Please avoid repeated freeze-thaw cycles.**

*The Buffer may vary depending on the Lot #. Please contact our technical support if you have specific requirements.

Amino Acid Sequence

APARSPSPST QPWEHVNAIQ EARRLLNLSR DTAAEMNETV EVISEMFDLQ EPTCLQTRLE LYKQGLRGLS
TKLKGPLTMM ASHYKQHCPP TPETSCATQI ITFESFKENL KDFLLVIPFD CWEPVQE

SDS-PAGE and Biological Activity of Human GM-CSF, cct-premium

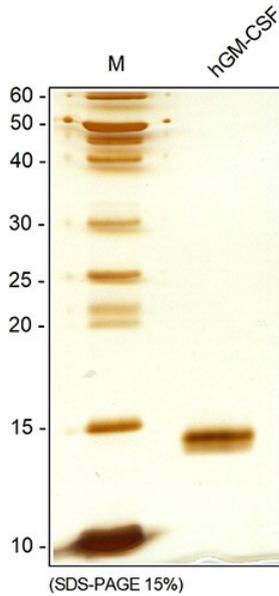


Fig. 1: SDS-PAGE analysis of recombinant Human GM-CSF, cct-premium. Sample was loaded in 15% SDS-polyacrylamide gel under reducing condition and stained with Silver stain.

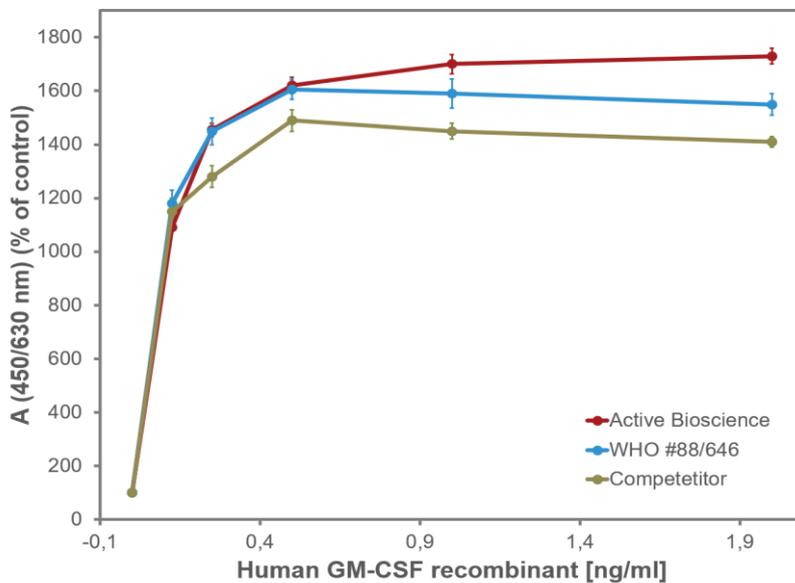


Fig. 2: Dose-dependent stimulation of cell proliferation in TF-1 cells by recombinant Human GM-CSF, cct-premium and the WHO standard 88/646. Values are the means (\pm SD) of triplicate determinations and expressed as percentage of control.

Usage: For research use only. Not for use in diagnostic or therapeutic procedures. Not for human use.

*The Buffer may vary depending on the Lot #. Please contact our technical support if you have specific requirements.

ORDERING
Tel.: +49 40 43208448-0
order@active-bioscience.de
www.active-bioscience.de

TECHNICAL SUPPORT
Tel.: +49 40 43208448-11
support@active-bioscience.de

Active Bioscience GmbH
Oberaltenallee 8
D-22081 Hamburg
HRB 98170 Amtsgericht Hamburg